
Introduction to ATAC-seq


Shamith Samarajiwa

University of Cambridge
CRUK Summer School in Bioinformatics
July 2018





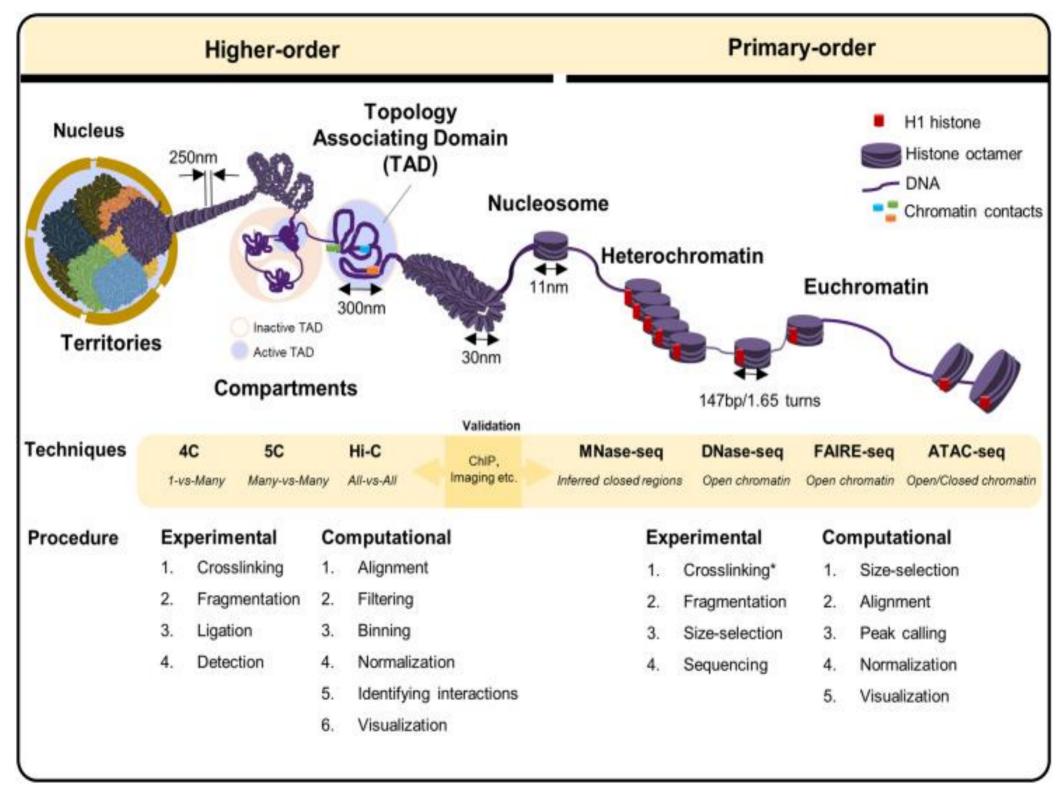

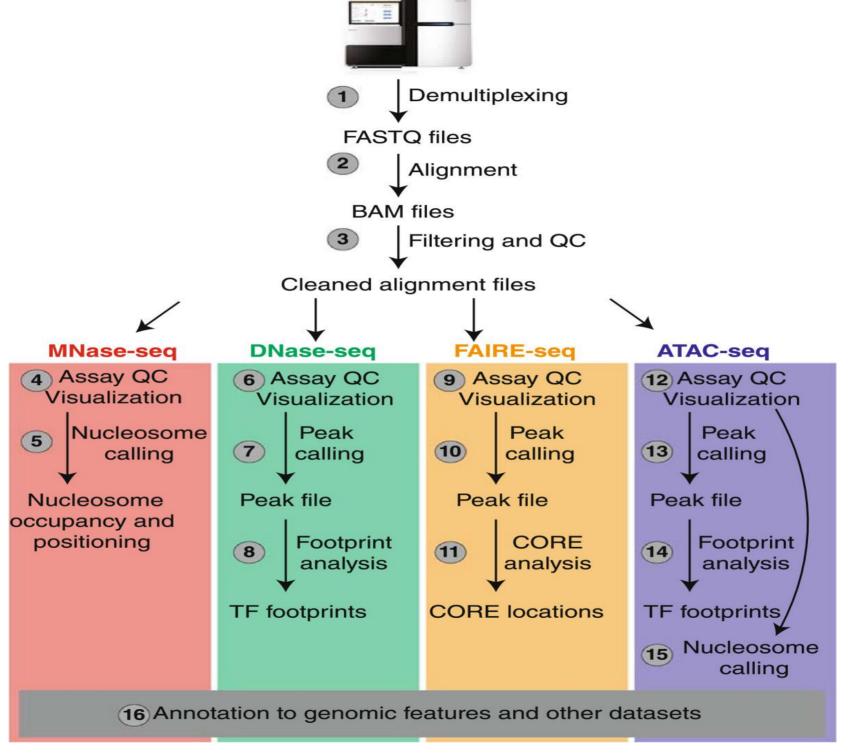

MENU 🗸


nature methods

Article | Published: 06 October 2013

Transposition of native chromatin for fast and sensitive epigenomic profiling of open chromatin, DNA-binding proteins and nucleosome position


Jason D Buenrostro, Paul G Giresi, Lisa C Zaba, Howard Y Chang [™] & William J Greenleaf [™]

Assay for Transposase Accessible Chromatin

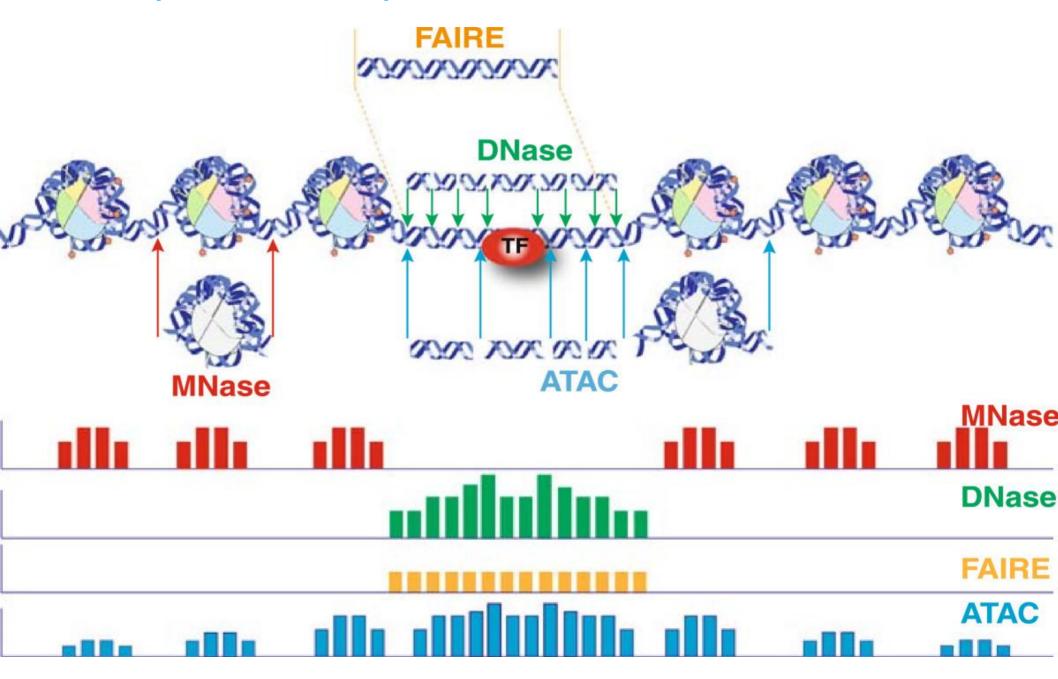


ATAC-seq


- Measure chromatin accessibility and nucleosome free regions
- Can also detect nucleosome packing, positioning and TF footprints
- Does not require sonication and phenol chloroform extractions, antibodies or sensitive enzymatic digestions that can introduce potential bias
- A hyperactive Tn5 transposase is used to fragment DNA and integrate into active regulatory regions
- During ATAC-seq, 500–50,000 unfixed nuclei are tagged in vitro with sequencing adapters by Tn5 transposase

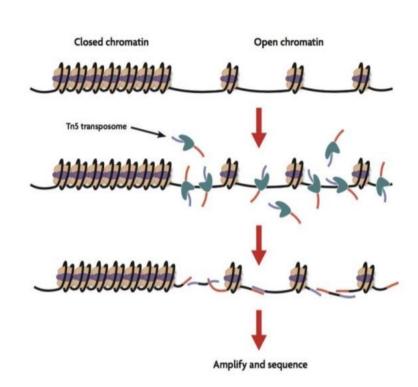


ATAC-seq

- Two-step protocol
 - Insertion of Tn5 transposase with adaptors
 - PCR amplification
- Needs ~500-50,000 cells
- Paired-end reads produce information about nucleosome positioning
- Insert size distribution of fragments has a periodicity of ~200 bp, suggesting that fragments are protected by multiples of nucleosomes
- Different fragmentation patterns can be associated with different functional states (eg. TSSs are more accessible than promoter flanking or transcribed regions)

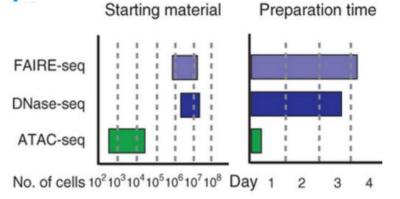


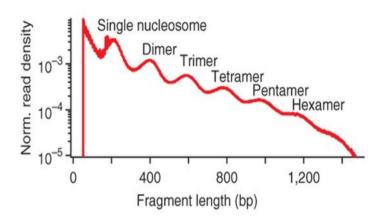

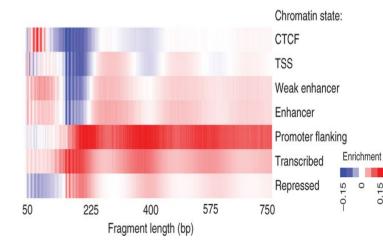

[1] JD Buenrostro et al, *Nature Methods*, 2013.

Workflow of ATAC-seq data processing

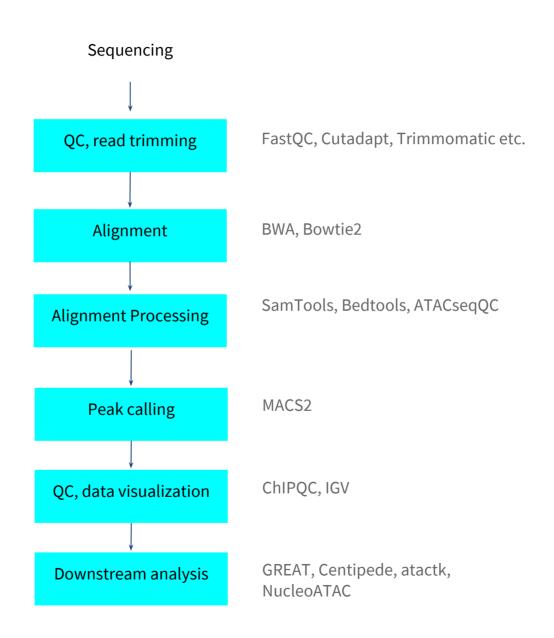

Working with ATAC-seq

- Remove mitochondrial reads
 - A large fraction of ATAC-seq reads map to the mitochondrial genome (up to 40-60%)
- Remove blacklisted regions before peak calling
 - hg19 or hg38 blacklisted regions can be obtained from ENCODE
- The signal (open chromatin where the transposase was inserted) is a mixture of feature types:
 - Short fragments signal from nucleosome free regions (NFRs) and open regions around DNA bound transcription factors. TSS rich.
 - Longer fragments open regions from around nucleosomes.
 Includes +1 and -1 nucleosome positions

Differences from ChIP-seq data processing

- Peak calling: Use fragment size for smoothing when calling peaks with MACS2
- Alternatively, to identify open regions following MACS parameters are used:
 MACS2 callpeak -t bamfile --nomodel --shift -100 --extsize 200 --format
 BAM -g hg38

For nucleosome occupancy shift and extension can centre the signal on nucleosomes (147 bp DNA is wrapped in a nucleosome)


MACS2 callpeak -t bamfile --nomodel --shift 37 --extsize 73 --format BAM -g hg38

Optionally, for paired end data

MACS2 callpeak -t bamfile --format BAMPE -g hg38

- --nomodel: don't build shifting model
- --shift: when this value is negative, ends will be moved toward 3'->5'
 direction
- --extsize: extend reads in 5'->3' direction to fix-sized fragments

ATAC-seq signal

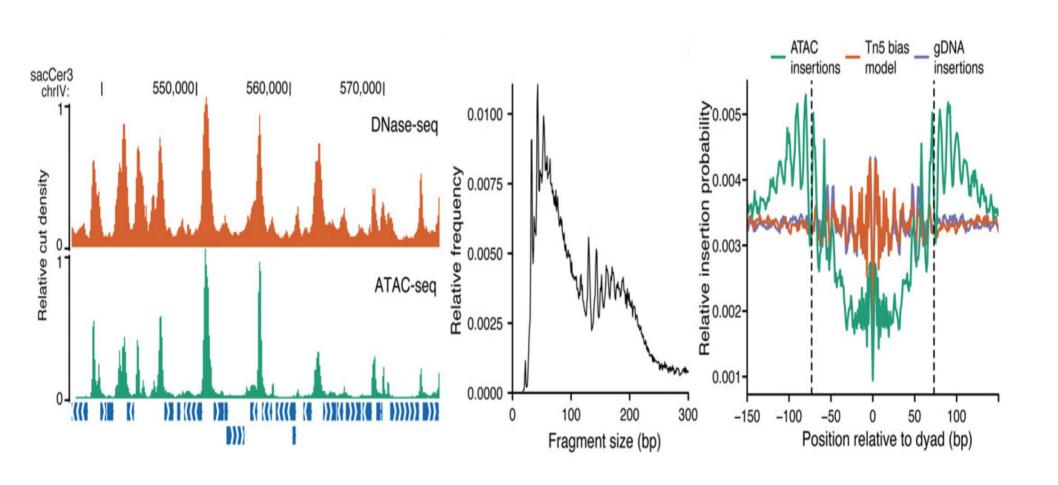

Normalization and Differential Accessibility of ATAC-seq

- Normalisation across samples might be needed
 - Efficiency of the ATAC-seq protocol in assaying open regions is affected by how many transposons get into the nuclei
 - One normalisation solution is to use signal from 'essential or housekeeping genes'. see: Sarah K. Denny et al, Cell, 2016.

- THOR is an HMM-based approach to detect and analyze differential peaks in two sets of ChIP-seq data from distinct biological conditions with replicates
 - Normalizes bam files (given a bed file of housekeeping genes)
 - Can compare two ATAC-seq datasets to perform differential openness analysis

Cutting site and TF Footprinting analysis

- ATAC-seq produces shorter fragments around smaller protected regions such as TF bound regions.
- Cut sites pileup signal is a good guide to TF binding sites
- Shifting reads +4/-5 nt depending on strand, should adjust for expected shift from Tn5 insertion.
- TF footprinting of these regions enables the detection of motifs for bound TFs
 - Centipede
 - msCentipede
 - o PIQ
 - Mocap